

Benchcap

Metallic fixed work bench. Equipped with two height adjusting jacks

Meuble de support fixe en métal. Equipé de deux vérins pour le réglage vertical

Tisch aus Metall, höhenverstellbar

Mueble de soporte fijo en metal. Equipado con dos elevadores para ajuste vertical

Mobile di supporto fisso in acciaio elettrozincato con piedini regolabili

Assembly Instructions Notice de Montage Instrucciones de Montaje Montage Beschreibung Istruzioni per il Montaggio













Necessary space and dimensions - according to the EN 14175-5 standard


Espace nécessaire et dimensions - selon la norme EN 14175-5

Abstände gemäß der Norm EN 14175-5


Espacio necesario y dimensiones – de acuerdo con la norma EN 14175-5

Spazio necessario e dimensioni - secondo la norma EN 14175-5

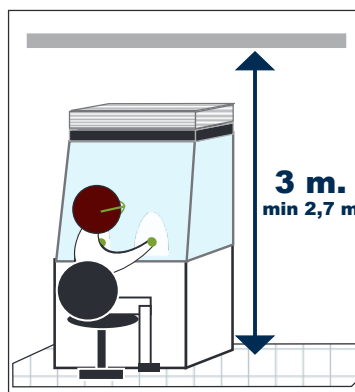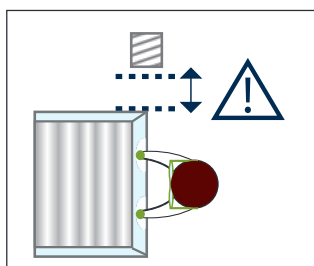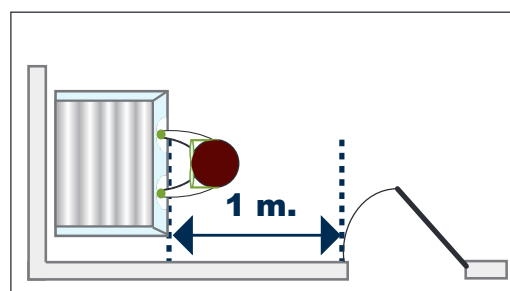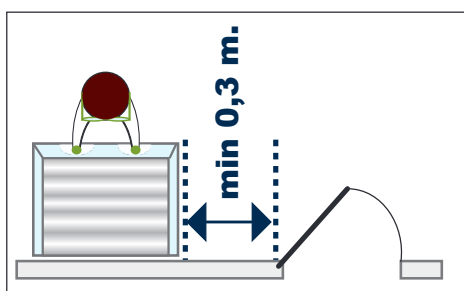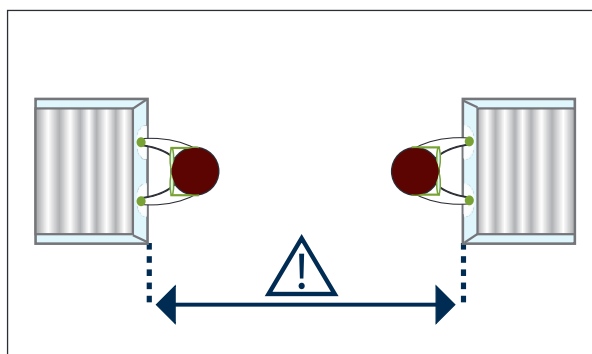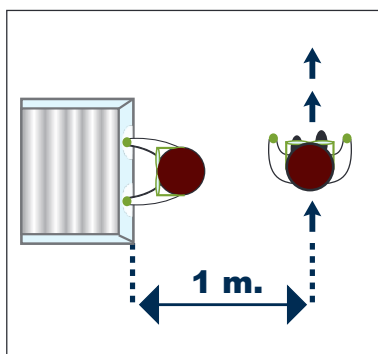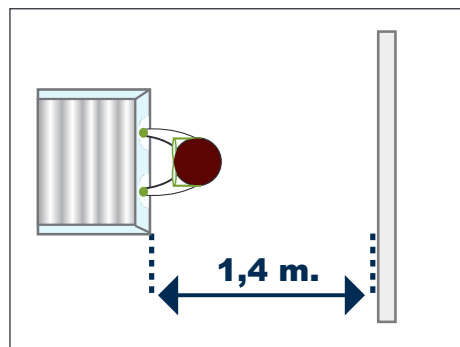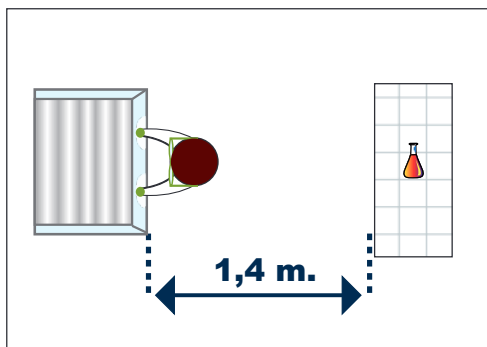


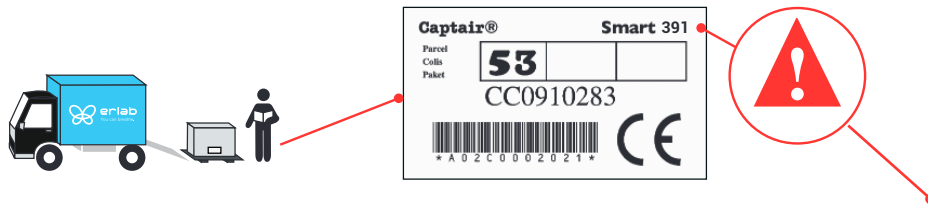


Example / Exemple / Beispiel / Ejemplo / Esempio : Benchcap 391 Smart

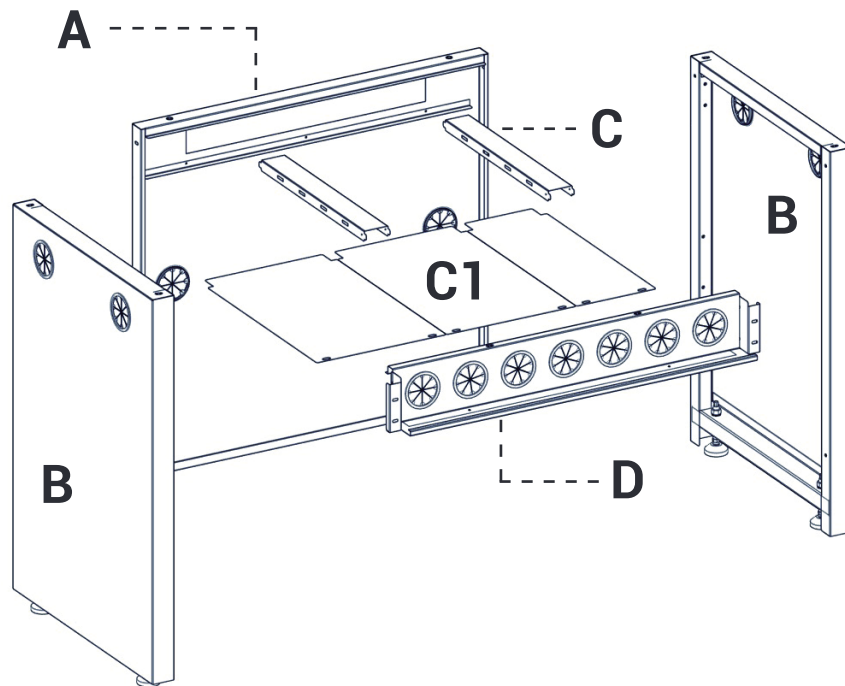


Benchcap 321

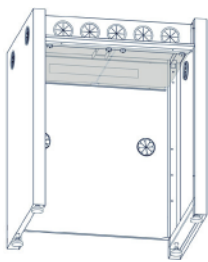


Benchcap 391/392

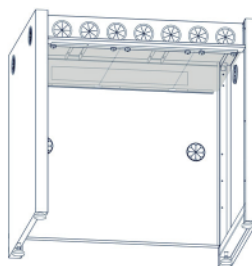


Benchcap 481/483



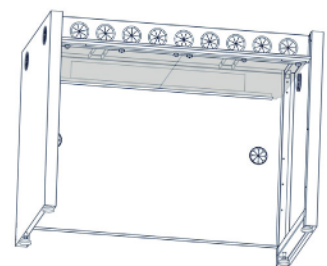


Benchcap 633

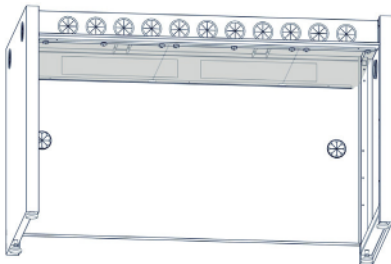


Benchcap 714

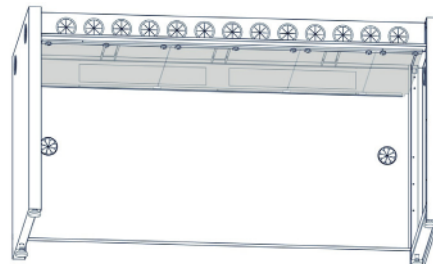

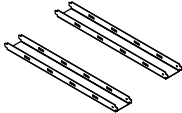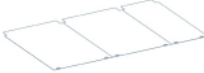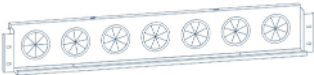
A	Back panel / Panneau arrière / Rückwand / Panel trasero / Pannello posteriore		x1	
B	Side panel / Panneau latéral / Seitenwand / Panel lateral / Pannello laterale		x2	
C	Crossbeam / Traverse / Querleiste / Travesaño / Traversa		321 x1 391/481/392/483/ 633 x2 714 x3	
C1	Board / Plateau / Plateau / Bandeja / vassoio		321/483 x2 391/481/392/633 x3 714 x4	
D	Utilities bar / Barre utilitaires		x1	
E	Nylon screws / Vis nylon / Nylonschraube / Tornillo nailon / Viti nylon		x2 + 321 x4 391/392/481/483/633 x8 714 x10	PIDB08547
F	Hex Head Screw / Vis tête hexagonale / Sechskantschraube / Tornillos de cabeza hexagonal / Viti con testa esagonale		x10	PIDB0X037
G	Rings / Rondelles / Scheibe / rondelle / arandelas		x10	PIDB0079
H	Stud / goujon / Gestüt / vite tassello / tornillo pasador		x2	PIDB0070
I	Optional shelf / Option étagère / Optionales Regal / Con estantería / Opzione mensola		x1	

1

A	x2
B	x1
F	x6
G	x6

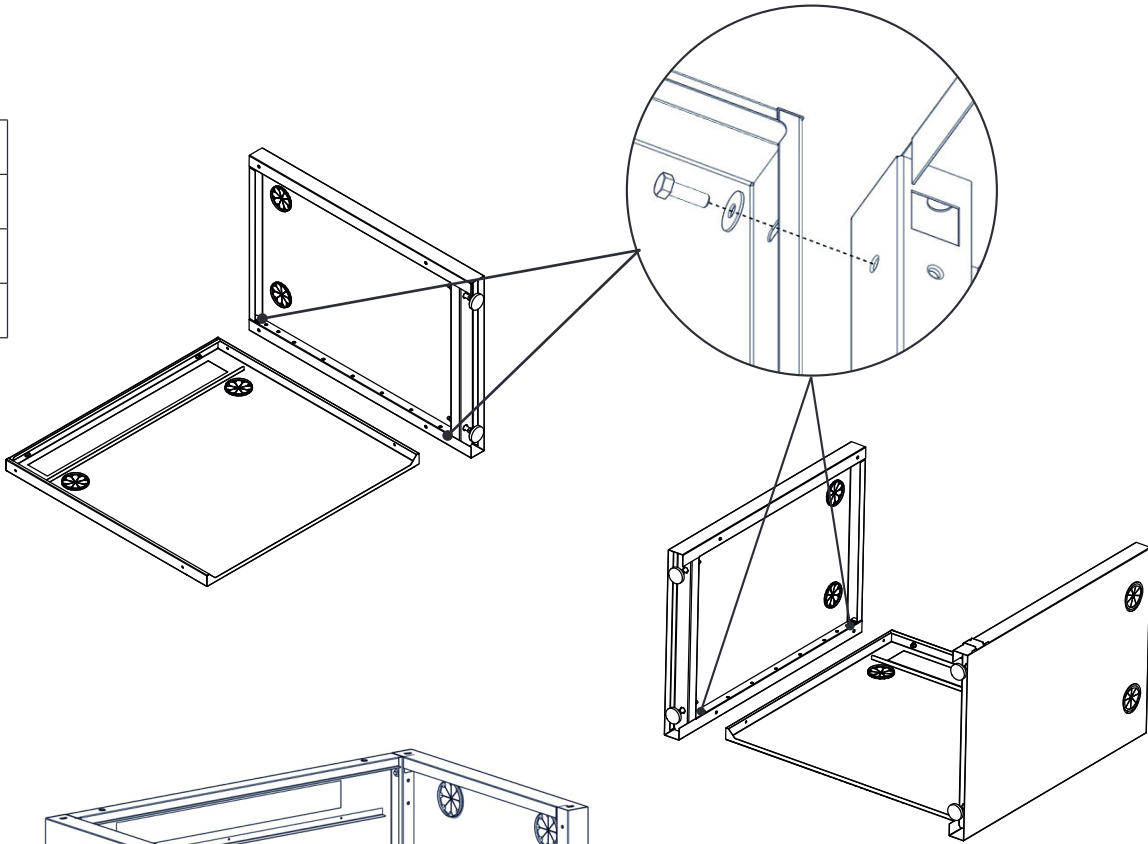


2

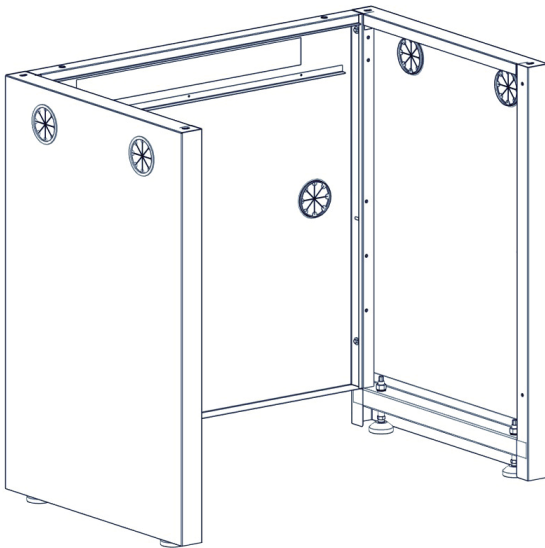


3

D	x1
F	x4
G	x4

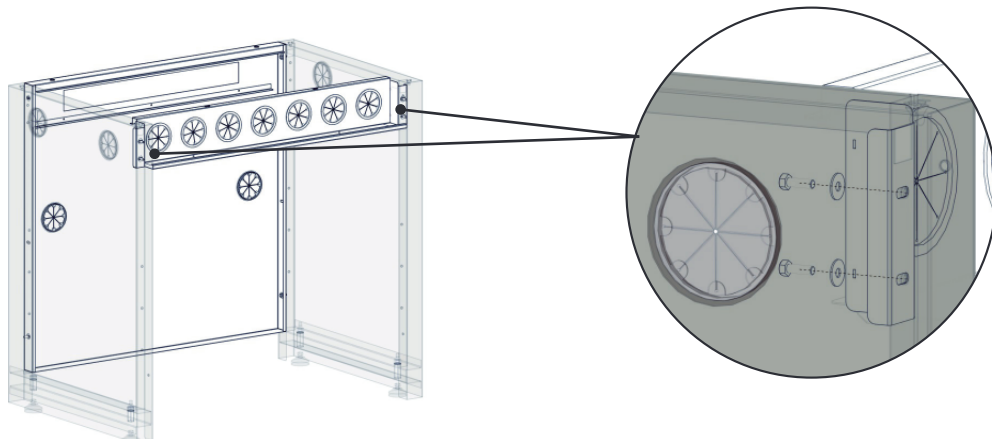


4

C1	321/483	x2
	391/481/392/633	x3
	714	x4
E	321	x4
	391/481/392/483/733	x6
	714	x8

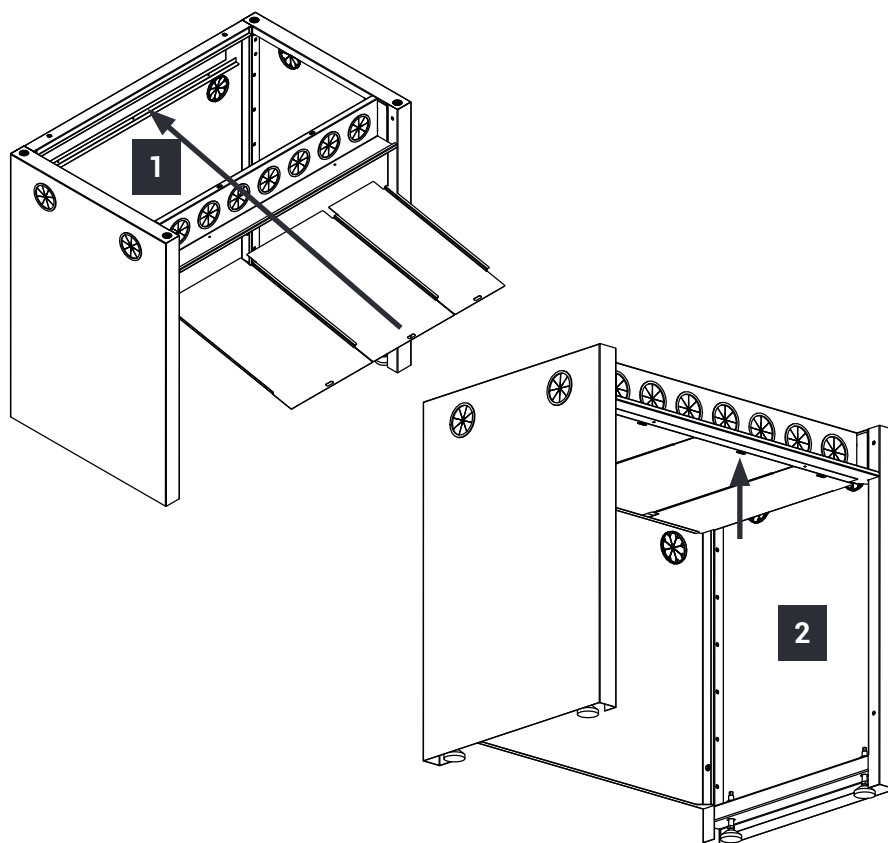


5

C	321	x1
	391/481/392/483/633	x2
	714	x3

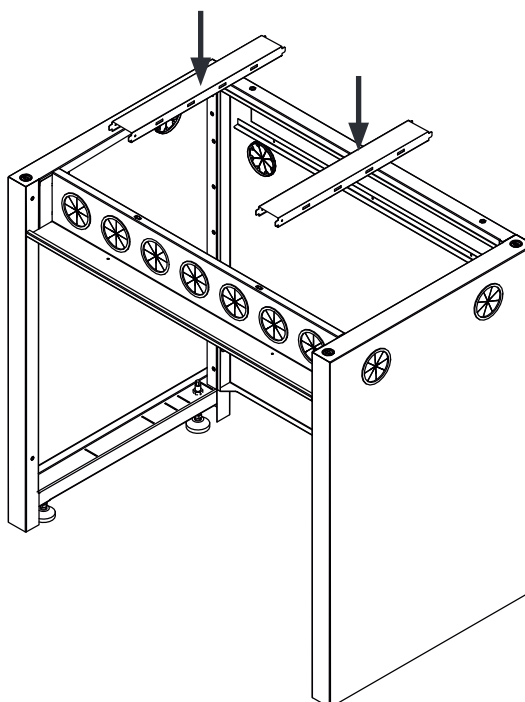

6



E x2

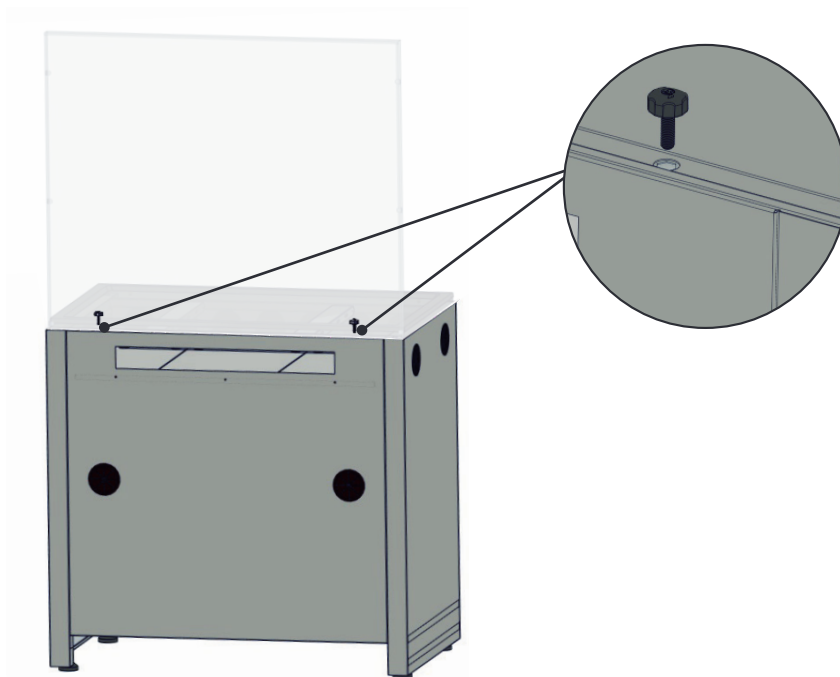

7

H x2

Work surfaces - Plan de travail - Piano di lavoro - Superficies de trabajo - Arbeitsplatten

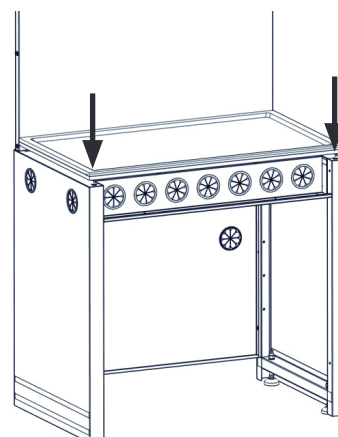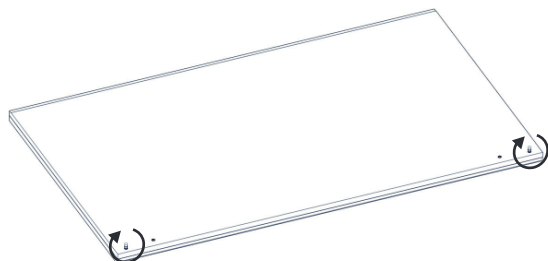

8

Optional shelf - Option étagère - Optionales Regal - Con estantería - Opzione mensola

I x1

